

# Chromation Spectrometer: Interface Information for Preliminary Datasheet

Mike Gazes

December 3rd, 2018

## Spectrometer Sequencing Diagrams

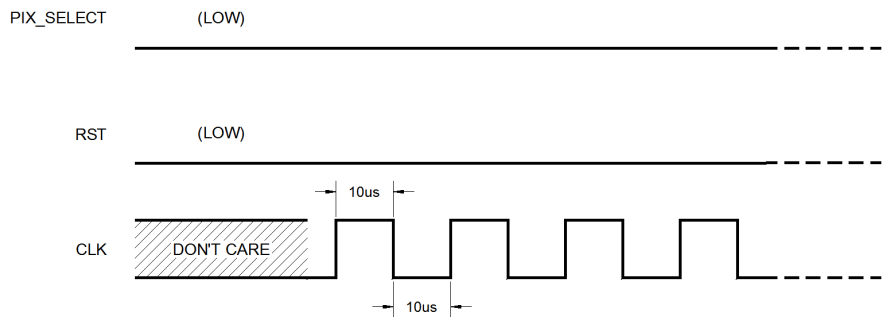


Figure 1: Spectrometer Initialization Sequence

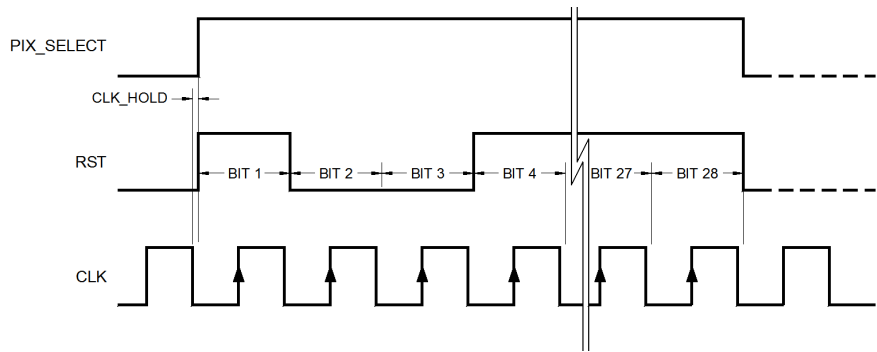


Figure 2: Spectrometer Configuration Sequence

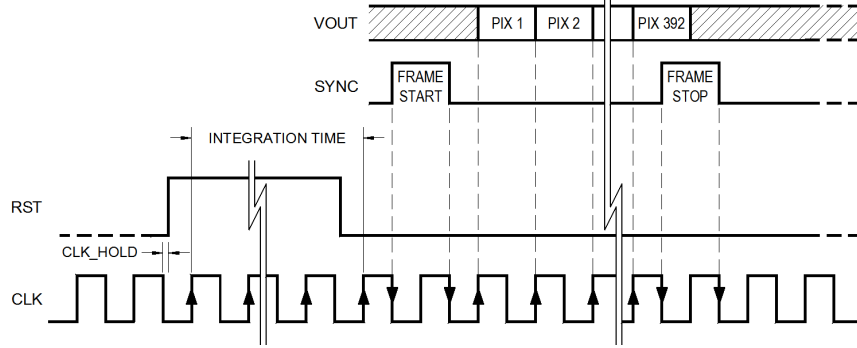


Figure 3: Spectrum Acquisition Sequence

## Recommended Operating Conditions

Table 1: Recommended operating conditions

Symbol	Parameter	MIN	TYP	MAX	Units
$V_{DD}$	Supply Voltage	2.8	3.0	3.3	V
–	Power Consumption while active	–	5.0	10	mW
–	Power Consumption in low-power mode	–	30	–	uW
$V_{in,HIGH}$	Input logic-level <i>HIGH</i>	$V_{DD}-0.7$	–	$V_{DD}$	V
$V_{in,LOW}$	Input logic-level <i>LOW</i>	–	–	0.7	V
$f_{CLK}$	Clock Frequency	15	50	200	kHz
$CLK_{HOLD}$	Clock Hold-time	–	10	–	ns
$CLK_{SETUP}$	Clock Setup-time	–	10	–	ns

Pins CLK, RST, and PIX\_SELECT are logic-level inputs. The spectrometer samples inputs on the **rising edges** of CLK.

$CLK_{HOLD}$  is the minimum time **after a falling edge of CLK** before RST is allowed to change value.  $CLK_{SETUP}$  is the minimum time RST must be stable **before a rising edge of CLK** to guarantee its value is sampled on the clock rising edge.

## Electrical Characteristics

Table 2: Electrical characteristics relevant to analog signal conditioning and ADC component selection

Parameter	Condition	MIN	TYP	MAX	Units
$V_{OUT}$ output impedance	–	–	10	–	kOhms
$V_{OUT}$ settling time	– 2	–	1	–	us
$V_{OUT}$ maximum swing	2.5x gain	–	$V_{DD}-0.3$	–	V
$V_{OUT}$ at dark	no light	0.60	0.84	1.1	V
Full well	pixels: 312.5um tall, 15.6um pitch	–	3.0e5	–	electrons
Linearity error per pixel	1x gain, $V_{OUT}=5\%-70\%$ full-well	0.5	1	5	%error
Conversion efficiency	pixels: 312.5um tall	–	6.5	–	uV/electron
Image lag	–	0.1	0.3	3.0	% $V_{SAT}$

is 0% error, indicating that an increase in the amount of photons incident on the pixel active area causes a linearly proportional increase in the dark-corrected output voltage.

*Full well* is the number of electrons a pixel stores at the saturation voltage. The saturation voltage,  $V_{SAT}$ , is the product of the *full well* and the *conversion efficiency*. When configured for 15.6 $\mu$ m-pitch, therefore,  $V_{SAT}$  is 1.95V and 70% of  $V_{SAT}$  is 1.365V.

Chromatation recommends dark-correcting the analog output before input to the ADC and selecting an ADC voltage reference between 70% and 100% of  $V_{SAT}$ . Of course,  $V_{SAT}$  and 70% of  $V_{SAT}$  do not correspond to standard voltage reference values. Voltage reference selection depends on application-specific design considerations.

## Recommended ADC values

Chromatation uses a 16-bit SAR ADC and 2.048V reference in the Chromatation Spectrometer Evaluation Kit. The voltage reference value is based on the Chromatation recommendations for linear pixel array configuration and analog dark-correction. See the reference design for the specific part numbers used in the Chromatation Spectrometer Evaluation Kit.

Table 3: Recommended ADC configuration for the spectrometer interface

Parameter	MIN	TYP	MAX	Units
Resolution	10	16	–	bits
Voltage reference	1.25	1.8	2.048	V

## Recommended Clock Filter

The rising edge of the CLK signal couples into the VOUT signal. Filter this transient with a simple passive RC filter with  $f_{3dB} \geq 1/2f_{CLK}$  or higher. This filter is not essential so it is safe to eliminate on ultra-low component-count designs. The table below shows recommended values.

Table 4: Recommended spectrometer VOUT passive RC clock filter

Symbol	Parameter	Condition	MIN	TYP	MAX	Units
$f_{3dB}$	Clock filter cutoff frequency	–	$1/2f_{CLK}$	–	–	kHz
$R_{filt}$	Clock filter R	$f_{CLK}=50kHz$	–	10	–	kOhms
$C_{filt}$	Clock filter C	$f_{CLK}=50kHz$	–	680	–	pF

Symbol	Parameter	Condition	MIN	TYP	MAX	Units
--------	-----------	-----------	-----	-----	-----	-------

## Spectrometer Configuration During Wavelength Calibration

Table 5: Spectrometer configuration used during wavelength calibration

Parameter	MIN	TYP	MAX	Units
Integration time	1	10	1000	ms
Internal analog gain	–	1x	–	–
Pixel height	–	312.5	–	um
Pixel pitch	–	15.6	–	um

## Spectrometer Interface

The following describes how to interface the Chromation spectrometer. Since the Chromation spectrometer is simply optical components mounted around the LIS-770i, *the spectrometer interface is the LIS-770i interface*. The interface description here is sufficient to operate the spectrometer but for simplicity some information about the LIS-770i is excluded. See the LIS-770i datasheet for additional operating conditions, waveforms, and details about other functionality not covered here, in particular the **power-down** mode.

External hardware is required to drive the spectrometer to achieve its higher-level function: *acquire a spectrum*. Chromation refers to this external hardware as the **spectrometer digital interface**.

### Spectrometer Interface Block Diagram

The block diagram below shows the spectrometer signals referred to in the following sequences to initialize, configure, and acquire a spectrum. The block diagram also shows how the **spectrometer digital interface** converts the spectrometer to a *SPI slave* device.

The purpose of the **spectrometer digital interface** is to hide the operational details of the LIS-770i and ADC, and instead provide the measurement system with a standard interface. Chromation recommends presenting the spectrometer interface to upstream devices as a *SPI slave*. A *SPI master* in the upstream

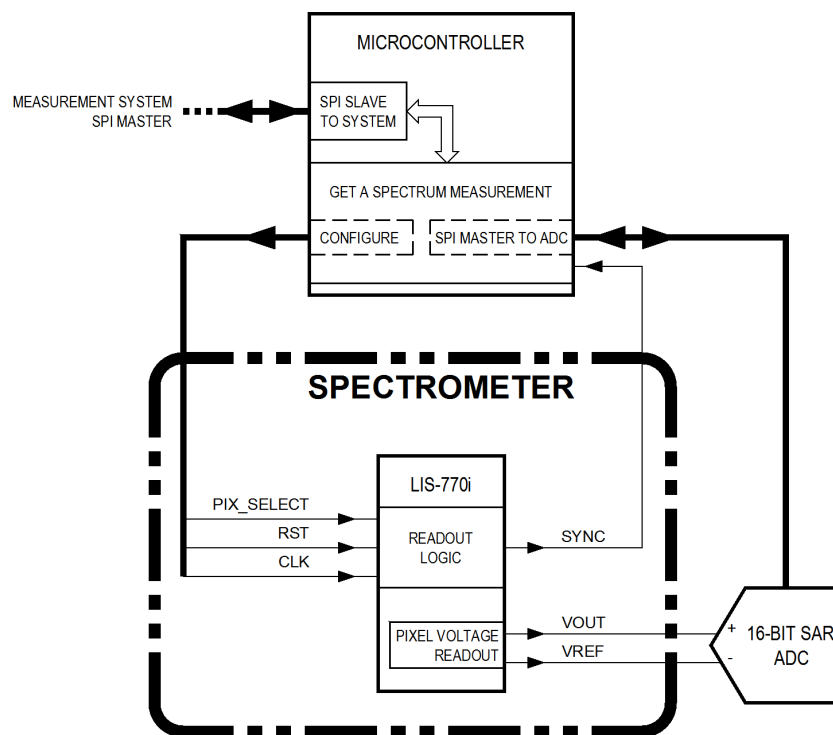


Figure 4: Spectrometer Interface Block Diagram

measurement system executes requests from the system host application by writing high-level commands recognized by the *SPI slave*.

## Initialize the spectrometer after power-up

The following sequence shows the recommended input logic-levels to apply after power-on. Start the clock signal on CLK after PIX\_SELECT and RST are driven *LOW*.

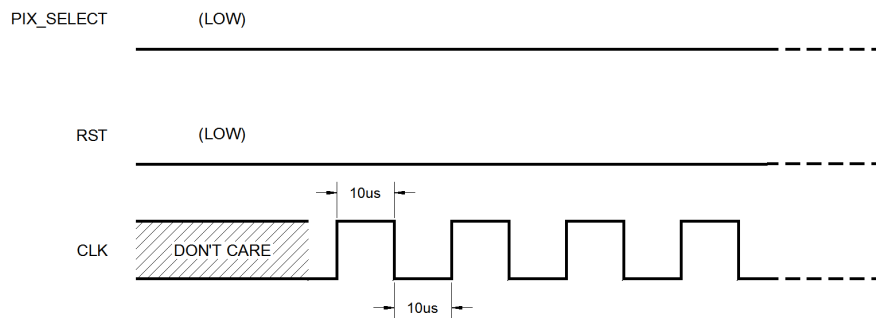


Figure 5: Spectrometer Initialization Sequence

- drive PIX\_SELECT *LOW*
  - PIX\_SELECT idles *LOW*
- drive RST *LOW*
  - RST idles *LOW*
- drive CLK with a 50kHz PWM, 50% duty cycle
  - the choice of clock speed is determined by the ADC conversion settling time and the speed of communication with the ADC
  - 50kHz is slow enough for a 16-bit SAR ADC to convert pixel voltages and to read out the converted digital value

## Configure the spectrometer

The spectrometer's linear pixel array has configurable parameters:

- analog gain: 1x, 2.5x, 4x, 5x
- pixel height: 62.5um, 125um, 187.5um, 250um, 312.5um
- pixel pitch: 7.8um (array is 784 pixels), 15.6um (array is 392 pixels)

The linear pixel array does *not* have a known default configuration on power-up. At a minimum, these parameters must be hard-coded in firmware and the linear pixel array must be programmed each time the device is powered-up.

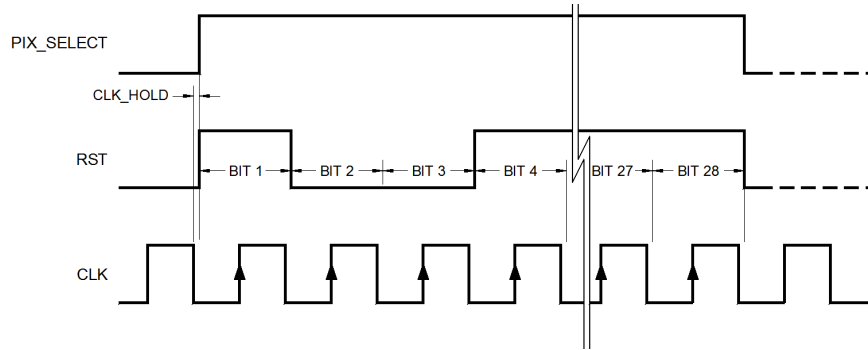


Figure 6: Spectrometer Configuration Sequence

Making these parameters configurable by the upstream measurement system may be beneficial in some context, particularly the analog gain. The spectrometer performance specifications only apply for the configuration intended by Chromation. This configuration matches the optical design of the Chromation spectrometer.

This is the intended configuration:

- analog gain: 1x
- pixel height: 312.5um
- pixel pitch: 15.6um (the array contains 392 pixels)

The following programming sequence configures the parameters of the linear pixel array parameters to match the intended configuration. The programming sequence consists of shifting in 28 bits. Bits are shifted in on the rising edge of CLK. The sequence starts with bit 1 and ends with bit 28.

- start the programming sequence:
  - drive PIX\_SELECT *HIGH* on a falling edge of CLK
- bit 1 *HIGH* programs the pixel pitch for 15.6um
  - drive RST *HIGH*
  - wait for the next falling edge of CLK
- bits 2 and 3 *LOW* set the analog gain to 1x
  - drive RST *LOW*
  - wait for the next falling edge of CLK
  - drive RST *LOW*
  - wait for the next falling edge of CLK
- the next 25 bits are *HIGH* to select the 312.5um pixel height
  - drive RST *HIGH*
  - wait for the next falling edge of CLK
  - repeat 24 more times
- stop the programming sequence:
  - drive RST *LOW*

- drive `PIX_SELECT LOW`

## Acquire a spectrum

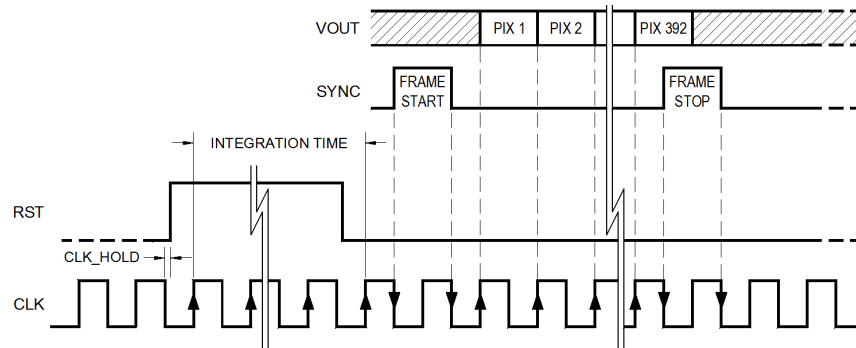


Figure 7: Spectrum Acquisition Sequence

The linear pixel array has internal switching logic to execute the spectrometer **exposure** and **readout** sequence, but like most other linear pixel arrays, the LIS-770i requires external logic-level signals to drive this internal switching. The LIS-770i also requires an external ADC to convert its analog output.

## Expose the spectrometer pixels

- drive `RST HIGH` to start exposure
  - drive `RST HIGH` after a falling edge of `CLK`
  - pixel exposure starts when the `HIGH` on `RST` is sampled on the rising edge of `CLK`
- wait the desired integration time by counting falling edges of `CLK`
  - for a 50kHz `CLK`, falling edges represent 20us ticks
  - for example, waiting 500 ticks is a 10ms integration time
- drive `RST LOW` to stop exposure
  - again, the `LOW` is sampled on the next rising edge of `CLK`

## Immediately readout the pixel voltages

- wait for a pulse on `SYNC`:
  - after `RST` is driven `LOW`, spectrometer output pin `SYNC` goes `HIGH` on the next falling edge of `CLK`
  - `SYNC` goes `LOW` on the next falling edge of `CLK` (`SYNC` is `HIGH` for one clock period)
- pixel readout begins on the next rising `CLK` edge after `SYNC` goes `LOW`



- start each ADC conversion on each rising edge of **CLK**
- finish each ADC conversion before the falling edge of **CLK**
  - the falling edge of **CLK** may couple into the pixel output voltage, so it is preferable to finish the ADC conversion before the **CLK** falling edge
  - and it is good practice in general to avoid transients during an ADC conversion
  - this is easily achieved using a 50kHz clock
  - in fact, using *Linear Technology #LTC1864L* as the ADC, both the conversion and readout of the converted samples fits within the **CLK** high time
- on the last pixel readout, **SYNC** pulses *HIGH* again for one clock period, starting on the falling edge of **CLK**

## Reference Design

This is the spectrometer interface in the Chromation Spectrometer Evaluation Kit.

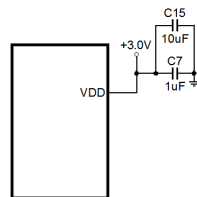


Figure 8: Reference Design